

RADINCON X-RAY Est 1979

Unit 19, 176 South Creek Rd, DEE WHY, NSW 2099
PO Box 5142, WHEELER HEIGHTS, NSW 2099

PH: 02 9971 7370
FAX: 02 9972 1734

ProImage - FILM ID PRINTER 9v

Light weight, Exposure adjustable, 9volt Battery. Small and compact in the darkroom. Comes with full instruction sheet

XA112 ~~\$270.00~~ + gst

SPECIAL
\$240.00+gst




FEBRUARY 2009

Daito XP 400 Auto Film Processor

No plumbing required • Easy to maintain • Automatic chemical replenishment • 3 minute processing
• Simple installation
• Perfect compact low volume unit

FINANCE AVAILABLE

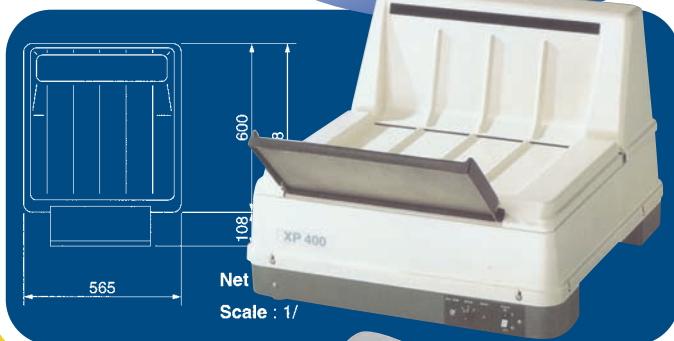

ProImage Cassettes & Screens



Ultra 60HF

High Frequency Portable X-ray Unit
3.2kW and only 12kgs.

Built especially for the equine vet looking for higher output this 3.2kW unit weighs in at only 12kgs • much the same weight as many other 2.0kW offerings. This is a powerful and compact option when it comes to the needs of larger in field work. The Ultra 60HF comes digital ready with a built in DR Interface.

Call for our combination package specials for this unit with the RAD-X mobile stands. This unit is set to become a class leader for this range of x-ray unit.



Phone: 1300 721 734


sales@radincon.com.au - www.radincon.com.au

QUALITY PRODUCT - RELIABLE SERVICE - HONEST ADVICE

A Division of RADINCON PTY LTD ABN 64 001 760 587

RAD-X HF PREMIER VET SYSTEM

The RAD-X HF Premier Vet Systems are the flexible choice when selecting a new x-ray suite for your hospital. We can tailor a package to suit your individual requirements and desires rather than an off the shelf system that may not be an ideal solution for you. The table tubestand combination provides the ultimate in positioning possibilities with movement in all practical directions. This is the best European manufactured offering we could find and surpasses the quality of all other options in the market. The RAD-X i and V Series High Frequency X-ray Generator choices are from 20 to 50kilowatt with the most popular veterinary being the 30kW and a matching x-ray tube with focal spots and performance to suit. Full Anatomical Programming is standard and in manual mode you have a choice of 2 or 3 point programming. This means you can select kV and mA or have more control of the image using kV, mA and time. These generators carry a 5 year manufacturers parts warranty. The RAD-X HF Premier Vet System is built to a performance, service and quality standard to keep you happy now and into the future.

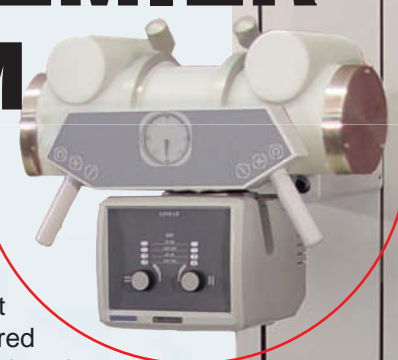

VTT Integrated Table/Tube Stand

Vet APR Console



FINANCE AVAILABLE


SE-30 RAD-X GEN STORED ENERGY

30kW High Frequency X-ray Generator

The SE-30 will provide all the x-ray energy needed for your Veterinary application. Plug it in to a 240V AC 15Amp outlet and you are ready to calibrate and make x-rays.

Utilizing the latest in high density electrolytic capacitor storage techniques, the SE-30 has no batteries to slowly loose their charge and require expensive replacement. Rated at more than 10,000 hours of continuous usage, these space age capacitors are at the heart of the stored energy solution provided by the SE-30.

No Batteries
30kW from 240V
15Amp supply
"The no compromise solution to the power problem."

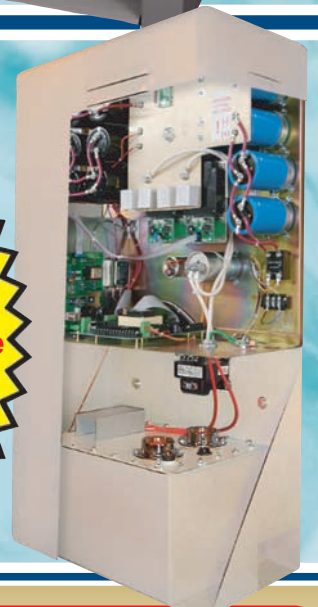

RADINCON X-RAY Est 1979

FEBRUARY 2009